Product Description

BIHOCL™ PureCleanse™ is a safe and effective irrigation, moisturizing, and debridement solution formulated to actively cleanse wounds and prepare the wound bed for healing. With no clinical contraindications for use, it can be used to prepare the wound bed for healing in a safe, effective and natural way.

Indications

BIHOCL™ PureCleanse™ is intended for cleansing, irrigating, moistening, debridement, and removal of foreign material including microorganisms and debris from exudating and/or dirty wounds, acute and chronic dermal lesions such as stage I-IV pressure ulcers, stasis ulcers, diabetic ulcers, post-surgical wounds, first and second degree burns, abrasions, minor irritations of the skin, diabetic foot ulcers, ingrown toenails, grafted and donor sites, and exit sites. It is also intended for moistening and lubricating absorbent wound dressings. The product is safe for use on adult and adolescent patients.

Ingredients

Water, sodium chloride, and 0.033% hypochlorous acid (330 PPM) as a preservative to inhibit growth of microorganisms within the solution.

Storage

Store the product in a room temperature environment and away from direct sunlight, heat, or extreme cold. Do not freeze. For adequate efficacy, do not use the solution past its expiration date (printed on the bottle itself).

Disposal

Dispose of remaining BIHOCL™ PureCleanse™ solution. Place bottle and sprayer in the regular waste disposal or biohazard waste disposal (if applicable). All components of the product packaging are recyclable. Please recycle wherever possible.

Instructions

To prepare for use, remove the sprayer from the protective pouch, remove the safety seal from the bottle, and firmly screw the sprayer onto the bottle.

1. Using a firm, consistent spray, direct the stream of solution onto the skin or wound, ensuring complete coverage.
2. Wait at least one minute. No rinsing is required.
3. No special disposal considerations are required.


Contraindications

Do not use on patients with known sensitivity to hypochlorous acid (HOCl).

Cautions

Not for injection. Single use only. Discontinue use immediately if irritation, sensitization, or reaction occurs. Keep out of reach of children.

This product is non-pyrogenic and is not made with natural rubber latex.

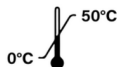Rx ONLY

MD

Health Canada Class II Medical Device License (MDL) #108007

Questions, comments or concerns? Please contact us.

info@biomiq.health | 1-800-573-1295 | www.biomiq.health

BIOMIQ

Manufactured for Biomiq Inc. by Sterasure Inc.

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Fabriqué pour BiomIQ Inc. par Sterasure Inc.